

HPF[®]



PORTABLE FOLDING MASSAGE TABLE

USER MANUAL

RETAIN THIS MANUAL FOR FUTURE REFERENCE

PLEASE READ THIS MANUAL CAREFULLY BEFORE USE


TABLE OF CONTENTS

PRECAUTIONARY INSTRUCTIONS..... 1
PARTS SCHEMATIC 2
SETTING UP YOUR TABLE..... 3
CLOSING AND TRANSPORTING YOUR TABLE..... 5

PRECAUTIONARY INSTRUCTIONS

1. Maintain the table periodically. Make sure that all screws are secure and the cable has not deteriorated or frayed. Make sure that there are no cracks on the aluminium structure.
2. Do not overload the table beyond the recommended weight limit. This could result in serious injury to the user. This product is tested to hold a maximum working load of 250 kg.
3. Do not use alcohol-based cleaners when cleaning the vinyl. This could cause damage to the vinyl material. Use warm, soapy water and wipe dry.
4. Make sure that when using the table, all four legs are at the same level. Using the legs at different levels can cause the table to be unstable and can cause injury.
5. Do not allow children to assemble or disassemble the massage table unless under adult supervision.
6. Be careful when using sharp objects, as they can pierce the covering. Keep the table's covering away from being pierced by the sharp items.

PARTS SCHEMATIC

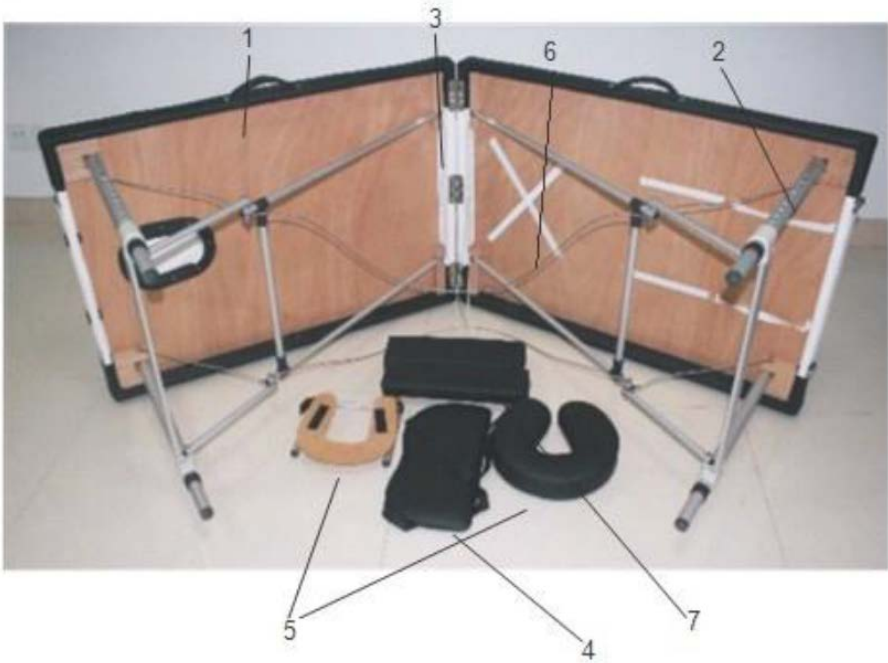

- 1 Handle
- 2 Copper button
- 3 Hinge
- 4 Hanged arm sling
- 5 Headrest
- 6 Cable
- 7 Side arm supports

SETTING UP YOUR TABLE

Steps are shown as the images below:



1



2



3



4



5



6

1. Release the buckle.
2. Lay the table on its side. With the four rubber feet on the floor, open up the table with the four legs extended as fully as possible.
3. Take out all accessories and then lift the table up until all of the legs are on the floor.
4. Put the headrest cradle into the corresponding outlets, position the cradle at an ideal angle for the face to sit comfortably. Lock in position by pressing the lever mechanism down (see figure below). Then put the headrest cushion on the cradle.



5. Remove the male clip attachments from the straps on the hanging arm rest. Re-attach them by inserting the strap through the underside of the clip's front slot, and down through the upper side of the back slot (see figure below). Then attach the hanged armrest or sling onto the headrest cradle by looping the straps through the bar and connecting the clips.



6. Insert the two side arm supports into the two holes at each side of the table.
7. The table is now ready for use.

The height of table is adjustable by pressing the copper buttons on the legs. All four legs should be done in the same way one by one and make sure that all of the legs are adjusted to the same height level.

CLOSING AND TRANSPORTING YOUR TABLE



1. Take out all accessories.
2. Lay the table on its side with the four rubber feet on the ground.
3. Put all accessories under the elastic straps under the top veneer.
4. Fold the table up and lock the buckle.
5. Put the massage table into the carry bag.
6. Carry the table like suitcase by grasping the handle.

WARNING!

If the frame/legs/table structure is damaged in anyway, you should not use this table. Using this table when damaged may cause harm or serious injury. Seek assistance from a qualified personnel immediately.



Some experts believe the incorrect or prolonged use of almost any product could cause serious injury or death. For information that may reduce your risk of serious injury or death consult the points below and additionally, the information available at www.datastreamserver.com/safety

- Consult all documentation, packaging and product labelling before use. Note that some products feature online documentation which should be printed and kept with the product.
- Check product for loose / broken / damaged / missing parts, wear or leaks (if applicable) before each use. Never use a product with loose / broken / damaged / missing parts, wear or leaks (if applicable).
- Products must be inspected and serviced (if applicable) by a qualified specialist every 6 months assuming average residential use by a person of average weight and strength, above average technical aptitude, on a property matching average metropolitan specification. Intended use outside these guidelines could indicate the product is not suitable for intended use or may require more regular inspection or servicing.
- Ensure all possible users of the product have completed an industry recognised training course before being given access to the product.

- The product has been supplied by a general merchandise retailer that may not be familiar with your specific application or your description of the application. Be sure to attain third party approval for your application from a qualified specialist before use regardless of prior assurances by the retailer or its representatives.
- This product is not intended for use where fail-safe operation is required. As with any product (take an automobile, aircraft, computer or ball point pen for example) there is always a small chance of a technical issue that needs to be repaired or may require replacement of the product or a part. If the possibility of such failure and the associated time it takes to rectify could in any situation inconvenience the user, business or employee or could financially affect the user, business or employee then the product is not suitable for your requirements. This product is not for use where incorrect operation or a failure of any kind, including but not limited to a condition requiring product return, replacement, service by a technician or replacement of parts could cause a financial loss, loss of employee time or an inconvenience requiring compensation.
- If this item has been purchased in error considering the points above simply contact the retailer directly for details of their returns policies if required.

